


### How to Operate:

- Plug power adapter plug into wall outlet and then plug power adapter into FIR Lumbar Support Pillow. Turn power switch on.
- When finished, turn power switch off. Disconnect power adapter from wall outlet and unplug pillow power adapter.

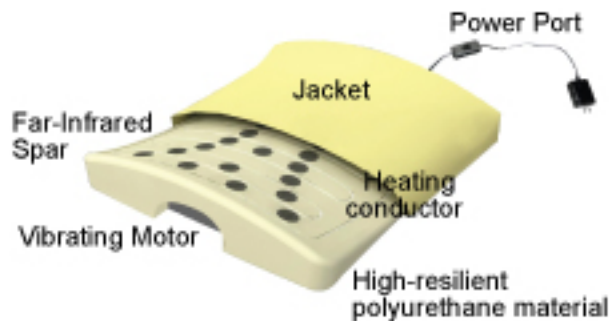

### Washing FIR Lumbar Support Pillow Jacket:

- Separate power adapter and FIR Lumbar Support Pillow.
- Remove plastic and metal nuts from FIR Lumbar Support Pillow power outlet turning counter clockwise using wrench.
- Pull zipper and fold inward to remove jacket.



### NOTES

- Do not fold or crease with use.
- Do not switch on power and use product when wet.
- Children should be under direct supervision.
- Do not use product if you are pregnant, wounded, have skin inflammation, skin diseases and/or cardiovascular diseases.
- Do not expose product to direct sunlight or high temperatures (no drying or ironing).
- Do not puncture product – keep away from sharp objects.
- Store product in a cool dry place.
- Do not place weight on top of unit.
- Note that cushion's jacket is washable, not the actual cushion.
- This product is not intended for those with physical/perceptual/mental impairments, or lack of knowledge or experience with those.

### Specifications

- Power Supply: DC 12V Input
- Power: <20W
- Parameter: Revolution 4000rpm
- Temperature 108-113°F
- Dimensions: 13" x 14" x 3" (L x W x H inches)
- Accessories: Auto DC Car Adaptor / Charger



*These statements have not been evaluated by the Food and Drug Administration, This product is not intended to diagnose, treat, cure or prevent any disease.*



HSIN TEN ENTERPRISE USA, INC. All rights reserved



## FIR Lumbar Cushion

Orthopedic Spinal Support Pillow  
w/ Massage and Far Infrared

**The FIR Lumbar Cushion** is an ergonomically designed lumbar support pillow made with high resilient polyurethane material and provides far infrared technology, vibrating massage and penetrating heat. The pillow conforms to the lower back to align posture, help prevent and relieve back pain and maintain kidney function.

**The FIR Lumbar Cushion** is equipped with Far Infrared strips on each side of the pillow. Lumbar support with heated massage and far infrared work together providing pain relief, corrected posture and optimized kidney function.

### Spinal Conditions that Cause Neck and Back Pain




Arthritis of the Spine - Bone Spurs - Bulging Disc - Degenerative Disc Disease - Herniated Disc - Pinched Nerve Facet Disease - Stenosis - Sciatica - Spinal Stenosis

### FEATURES:

#### Far Infrared Rays (FIR)

The FIR Lumbar Cushion utilizes far infrared technology to expand capillaries which assist in increased blood flow and circulation.

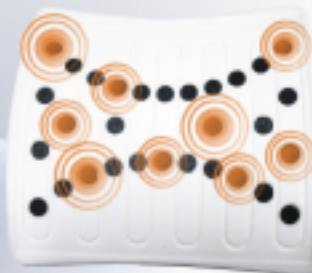

#### Heating Conductor

The FIR Lumbar Cushion features a heating conductor that offers soothing relief to sore muscles and joints and can reach a temperature of up to 113°F.




#### Vibrating Massage

The FIR Lumbar Cushion consists of two individual motors to help relieve lower back pain, improve circulation and proper function to the kidneys.



#### Ergonomic Design

The FIR Lumbar Cushion was specifically engineered to adequately support the muscles and skeletal structure of the lower back and supports good posture in individuals.

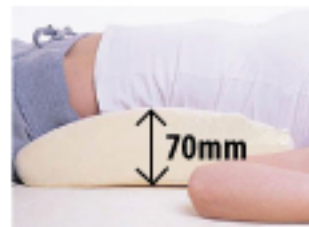

### A Combination of Excellence.

The FIR Lumbar Cushion is equipped with Far Infrared strips on each side of the pillow. Lumbar support with heated massage and far infrared work together to provide pain relief, correct posture and optimize kidney functioning.

### The FIR Lumbar Cushion is Suitable for...



#### Office Workers

Ideal for those who spend long hours behind a desk.



#### Women

Aids in the relief of pain and discomfort associated with menstrual cramps.

Ideal for post pregnant mothers. Provides soothing relief to the pelvic region.



#### Elderly

Relieves muscle pain caused by poor blood circulation and inactivity.



#### Students and Teenagers

Supports good posture and proper spinal alignment during adolescence.



#### Drivers

Provides an excellent source of relief while driving. Ideal for long commutes and road trips.



#### Athletes and Laborers

Helps relieve soreness and lower back pain after physical training or hard labor.